Breeze[™] Heating/Cooling Work Station

When combined with a circulator, the compact Breeze provides rapid heating/cooling and is ideal for applications requiring precise control by solution temperature

Designed as an add-on module for the Carousel 6 Plus and Tornado, Breeze creates a parallel process reactor for controlled heating and cooling. Breeze is ideal for applications that require fast and precise solution temperature control, such as crystallisation studies.


Features

- Thermofluid -85 °C to +235 °C providing a solution temperature of -30 °C to +165 °C.
- 135 mm ø top plate integrates with the Carousel 6 Plus, Tornado, and Heat-On.
- Breeze's small internal volume ensures a quick response to changes in thermofluid temperature.



Storm[™] Heating/Cooling Work Station

Combined with a suitable circulator, Storm can provide controlled steady state heating and cooling

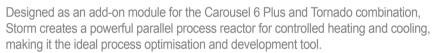

Features

- Thermofluid -85 °C to +235 °C providing a solution temperature -65 °C to +200 °C.
- 135 mm ø top plate integrates with the Carousel 6 Plus, Tornado and Heat-On.
- Unique internal design maximises heat transfer whilst the insulated outer case reduces heat loss and prevents contact with hot/cold thermofluid.





165°C

200°C

65°C

Breeze with 250 ml Heat-On, stand and overhead stirrer



Storm with Carousel 6 Plus, Tornado, overhead stirrer and PTFE insulating plate



Carousel 6 Plus locates on to the Storm without tools