



Cool-It is virtually unbreakable and designed to allow powerful stirring with all leading brands of stirrer including Heidolph, IKA, Eyela and Radleys.

## Easy pour spout and handle

Cool-It's unique non-drip spout and ergonomically designed handle makes the disposal of waste solvents much safer and easier than with a standard dewar or bowl; avoiding accidental spillages and creating a safer working environment.

## Large and small bowl options

Cool-It is available in two sizes of bowl. The small Cool-It bowl is suitable for flasks up to 400ml, and large Cool-It bowl for flasks up to 2 Litres. Both sizes have their own dedicated lid.

## Protects your stirrer and minimises spills

As well as insulating your reaction Cool-It minimises condensation on the outer surfaces, protecting your stirrer from frost and moisture ingress. It is also designed to fit securely onto the top plate minimising accidental spills and the risk of the bowl being knocked off the top plate.

## Cool-It fits on to all popular brands of stirring hotplate

Cool-It bowls are suitable for use with all stirring hotplates with a circular top plate up to 145mm in diameter. For stirring hotplates with 135mm ø top plates (e.g. IKA, Eyela, Radleys etc) an optional top plate adapter is available if required.

## **Cool-It accessories**

A useful range of accessories are available including, clamp, stand, stirring hotplate, digital thermometer, scoop, protective gloves, apron and face shield.





Large Cool-It bowl for flasks up to 2 litres





Large Cool-It bowl with lid, stirrer, stand and digital thermometer

# Cool-It<sup>™</sup> insulated bowls for round bottom flask up to 2 litres

The safest and most efficient way of cooling and stirring round bottom flasks to -78°C.



Small Cool-It bowl for flasks up to 400ml

Cool-It replaces fragile glass dewars, unstable plastic bowls and keeps your chemistry colder for longer.

The compact and virtually unbreakable Cool-It insulated bowls are designed to fit onto a standard stirring hotplate to cool and stir round bottom flasks, beakers, test tubes etc.

## Cool-It keeps it cooler for longer


The excellent design and insulation qualities of the Cool-It bowls keep your sample cooler for longer.

- Compared to a standard non-insulated plastic bowl Cool-It will keep your sample below -70°C for up to 5 times longer.
- Compared to a glass dewar Cool-It will keep your reactions below -70°C for twice as long.







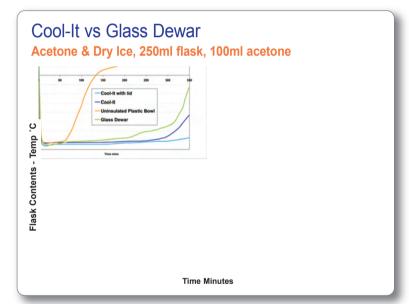

## Optional lid increases cooling time by 20%

Cool-It bowls come with unique two piece lid, which is easily fitted (and removed) once your flask is in place. Using the lid can keep your reaction cool for up to 20% longer, minimises condensation and ice formation on your flask (maintaining visibility of the contents) and prevents spitting from the cooling mixture.

### **Unbreakable**

Cool-It bowls are manufactured from a robust, chemically resistant HDPE casing encapsulating a high quality insulated foam core. The combination of tough composite materials not only provides excellent insulation but is also durable, and unlike fragile glass dewars, Cool-It is virtually unbreakable if dropped.



# Cool-It<sup>™</sup> Specifications and Ordering Information



#### **Cool-It Bowls**

#### Cool-It Small Bowls

 RR71011
 Cool-It Small Bowl + Free Lid

 RR71012
 Cool-It Small Bowl

 RR71015
 Cool-It Small Lid

RR61085 Adapter Plate for 135mm Stirring Hotplates

### Cool-It Large Bowls

RR71021 Cool-it Large Bowl + Free Lid
RR71022 Cool-it Large Bowl
RR71025 Cool-it Large Lid

RR61085 Adapter Plate for 135mm Stirring Hotplates

## **Cool-It Systems**

#### Cool-It Basic System 1

RR71500 Cool-It Basic System 1

1 x RR71011 Cool-It Small Bowl + Free Lid 1 x RR71021 Cool-It Large Bowl + Free Lid

1 x RR61085 Adapter Plate for 135mm Stirring Hotplates 1 x RR71200 Stirring Bar Evaluation Kit - for flasks (pk 10)

1 x RR99908 Dry Ice Scoop

### Cool-It Standard System 2

RR71501

Cool-It Standard System 2

1 x RR71011 Cool-It Small Bowl + Free Lid 1 x RR71021 Cool-It Large Bowl + Free Lid

1 x RR61085 Adapter Plate for 135mm Stirring Hotplates

1 x RR99905 Digital Thermometer (-250C to +400C) and 200mm Probe

1 x RR71100 Support Stand A-Shape Base 1 x RR71105 Support Rod 600x12.5mm 1 x RR71110 Retort Clamp to 85mm 1 x RR71115 Boss Head 16mm

1 x RR71200 Stirring Bar Evaluation Kit - for flasks (pk 10)

1 x RR99908 Dry Ice Scoop

## Cool-It Complete System 3

RR71502

Cool-It Complete System 3

1 x RR71011 Cool-It Small Bowl + Free Lid 1 x RR71021 Cool-It Large Bowl + Free Lid

1 x RR61085 Adapter Plate for 135mm Stirring Hotplates 1 x RR98072 Carousel Stirring Hotplate 230v UK Plug

1 x RR99905 Digital Thermometer (-250C to +400C) and 200mm Probe

1 x RR71100 Support Stand A-Shape Base 1 x RR71105 Support Rod 600x12.5mm 2 x RR71110 Retort Clamp to 85mm 2 x RR71115 Boss Head 16mm

1 x RR71200 Stirring Bar Evaluation Kit - for flasks (pk 10)

1 x RR99908 Dry Ice Scoop

1 x RR99909 Cold Temperature Apron 1060mm long, Waterproof

1 x RR99910 Protective Face Shield

1 x RR98024 Protective Cold Temperature Gloves

## Cooling Protection Kit

RR71505 Cooling Protection Kit

1 x RR99909 Cold Temperature Apron 1060mm long, Waterproof

1 x RR99910 Protective Face Shield

1 x RR98024 Protective Cold Temperature Gloves

## Clamp & Stand Kit

RR71510

Flask Stand & Clamp Kit

1 x RR71100 Support Stand Base A- Shape 1 x RR71105 Support Rod 600x12.5mm 1 x RR71110 Retort Clamp to 85mm 1 x RR71115 Boss Head 16mm

#### **Cool-It Accessories**

### Stirring Hotplates & Temperature Controller

RR98072 Carousel Stirring Hotplate 230v UK Plug
RR98072/EURO Carousel Stirring Hotplate 230v Euro Plug
RR98072/JAP Carousel Stirring Hotplate 110v Japanese Plug
RR98072/SWISS Carousel Stirring Hotplate 240v Swiss Plug
RR98072/USA Carousel Stirring Hotplate 115v US Plug
RR98073 Carousel Temperature Controller

#### PTFE Magnetic Stirring Bars

RR98075 Cross Shape Stirring Bar 10mm (for 5 to 10ml) pk 40

RR98091 Cross Shape Stirring Bar 16.5mm RE (for 25 to 50ml) pk 20

RR99064 Elliptical Stirring Bar 25mm RE (for 100 to 400ml) pk 10

RR95920 Oval Stirring Bar 40mm (for 500 to 1000ml) pk 5

RR95921 Oval Stirring Bar 50mm (for 2000ml) pk 5

RR71200 Stirring Bar Evaluation Kit - for flasks, pk 10

(contains 2 each of the above stirring bars)

#### Cool-It Accessories

RR99905 Digital Thermometer (-250C to +400C) and 200mm Probe

RR99906 Digital Thermometer (-250C to +400C)
RR99907 200mm Temperature Probe
RR99908 Dry Ice Scoop

RR99909 Cold Temperature Apron 1060mm long, Waterproof

RR99910 Protective Face Shield
RR98024 Protective Cold Temperature Gloves

## Support Stands & Clamps

 RR71100
 Support Stand A-Shape Base

 RR71105
 Support Rod 600 x 12.5mm

 RR71110
 Retort Clamp to 85mm

 RR71115
 Boss Head 16mm

Prices are in Pounds Sterling or Euros and exclude VAT and delivery. For pricing in other currencies please contact your local distributor

